

Hi-Flo MiniHEART® Circulaire® II

Continuous Nebulizer

Optimize
efficiency,
economy and
clinical effect.



Westmed

When Your Care Is Critical

Hi-Flo MiniHEART® Circulaire® II

Continuous Nebulizer

Realize the Best of Both Worlds.

Now you can achieve the efficiency and positive effects of continuous nebulization therapy (CNT) with the advantages of enhanced delivery, less waste and greater protection. It all comes together in the new Hi-Flo MiniHEART® Circulaire® II Continuous Nebulizer, exclusively from Westmed.


Start with the proven performance of the MiniHEART® continuous delivery nebulizer that eliminates the inconsistencies associated with intermittent aerosol delivery, while reducing valuable clinician time required to administer therapy. Now add the innovative Circulaire® II delivery system with its 550 mL reservoir bag that captures aerosolized medications during patient exhalation, to boost delivery upon subsequent inspirations, while reducing drug waste and minimizing room air entrainment. Plus, the variable expiratory resistor enables Positive Expiratory Pressure (PEP) to enhance drug deposition. This unique combination of features makes the new Hi-Flo MiniHEART® Circulaire® II Continuous Nebulizer the ideal choice in your ER.

Select the Lo-Flo MiniHEART® Continuous Nebulizer for medication delivery during mechanical ventilation and other situations that require no more than 2 L/min of additional flow rate.

Protect the Patient. Protect the Caregiver.

Advanced safety features of the new Hi-Flo MiniHEART® Circulaire® II Continuous Nebulizer include an integral expiratory filter that protects the caregiver from inhaling medications or patient exhalation. Also, the Circulaire's integral one-way valve virtually eliminates the risk of the patient exhaling droplets into the device, to protect both the nebulizer and the reservoir from contamination by exhaled organisms. In addition, the Comfort Soft Plus™ non-vented mask ensures a snug and comfortable fit throughout the therapy.

- Timely, accurate delivery of continuous aerosolized medication
- 30 mL MiniHEART® nebulizer allows up to 60 minutes of therapy per fill
- Closed system minimizes air entrainment & provides high FIO₂
- Enhanced medication delivery performance with less waste
- Unequalled patient and caregiver protection
- Shorter transition time from impending respiratory failure to stability
- Quicker return to patient baselines
- Less time in the ER
- Lower hospital admission rates
- Improved labor efficiency

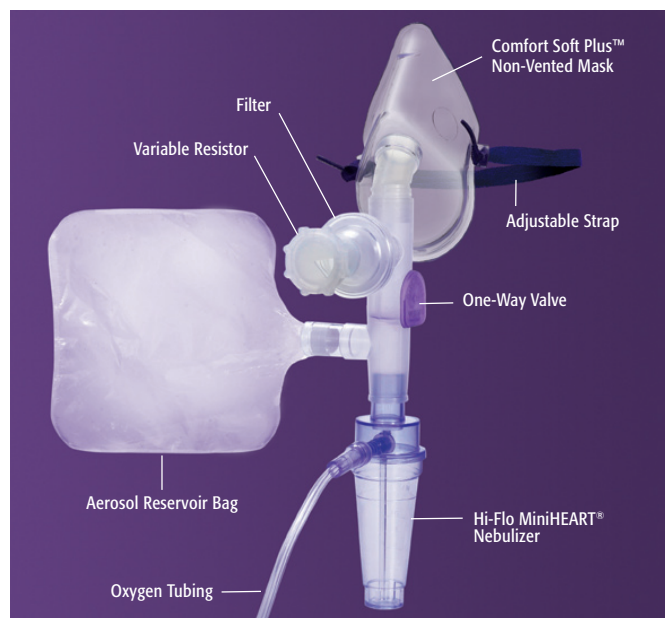

Part Number	Description	Case Quantity
100339	Hi-Flo MiniHEART® Circulaire® II Continuous Nebulizer with Adult Comfort Soft Plus™ Non-Vented Mask and 7' Tubing	25
100612	Hi-Flo MiniHEART® Nebulizer	25
100612A	Hi-Flo MiniHEART® Nebulizer with Adult Mask and 7' Tubing	25
100613A	Hi-Flo MiniHEART® Direct Connect with Pediatric Mask and 7' Tubing	25
100614A	Hi-Flo MiniHeart® with Pediatric Dragon Mask and 7' Tubing	25

Westmed

When Your Care Is Critical

Westmed, Inc.

5580 S. Nogales Highway
Tucson, AZ 85706
Toll Free: 800-975-7987
Phone: 520-294-7987
Fax: 520-294-6061
www.westmedinc.com

Distributed by:



Visit www.westmedinc.com
for a complete listing of
Westmed products.

US and International Patents Pending.
©2014 Westmed, Inc. All rights reserved. PN 75310, Rev .01